



Easy to apply anterior closure

12"



16"

KNAPP™ HINGED KNEE ORTHOSIS - Anterior Closure

- Constructed of Kuhl™ perforated neoprene for cool, comfortable compression
- 2 Lengths effectively handle multiple instability diagnoses
- 3 Hinge styles increase customization
- Hinges are fully covered for sports activity, yet easily accessible
- Elasticized band in thigh strap ensures comfort
- Ideal for pre/post-operative care of MCL/LCL and meniscal repairs
- Longer length can be used for mild to moderate ACL/PCL sprains & strains and as an adjunct to patella fracture



Thomas P. Knapp, M.D.

Developed in collaboration with Thomas P. Knapp M.D., Santa Monica Orthopaedic and Sports Medicine Group, Santa Monica, CA.

Order No.	Description / Options
5656	Kuhl™ Knapp™ -12" Length, Standard Hinges
5656 HH	Kuhl™ Knapp™ -12" Length, Heavy Hinges
5656 ROM	Kuhl™ Knapp™ -12" Length, ROM Hinges
5658	Kuhl™ Knapp™ -16" Length, Standard Hinges
5658 HH	Kuhl™ Knapp™ -16" Length, Heavy Hinges
5658 ROM	Kuhl™ Knapp™ -16" Length, ROM Hinges



Open popliteal eliminates bunching, improves range of motion and cools the back of the knee.

Hinges are offset in the center to reduce condyle irritation



3 Hinge Options:

1. Flip Top Pin Stop Range Of Motion Hinge
2. Standard Polycentric Hinge
3. Heavy Polycentric Hinge

ROM Hinge flexion stop placement:

In 0°, 10°, 30°, 50°, 70°, 90° and 110° settings.

ROM Hinge extension stop placement:

In 0°, 10°, 30°, 50°, 70°, 90° and 110° settings.


The ROM hinge will also allow you to set it in full extension. (It should be noted that due to neoprene's four-way stretching character, care should be taken to insure that your patient does not drive through your pre-set stop.)



Neoprene condyle pads provide protection and comfort and reduce chaffing.




Moisture Escapes

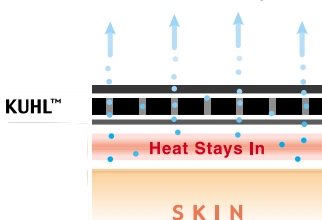

Kuhl™ Neoprene - 1/8" perforated neoprene laminated between 2 nylon layers.

The perforations allow moisture to escape, allow the skin to remain dry, while acting as an insulator of body heat.

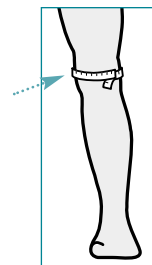

KUHL™ KNAPP™ KNEE SIZING:

Measure at mid-patella with leg fully extended.

- | | |
|---------------------|---------------------|
| XS: 10"-12" | Sm: 12"-14" |
| Med: 14"-16" | Lg: 16"-18" |
| XL: 18"-20" | XXL: 20"-22" |

Larger and custom sizes are available upon request.