


**Warnings and Instructions: Review carefully, proper application is required**


**⚠ Warning:** This device will not prevent or eliminate risk of injury. Do Not Overtighten. If swelling, pain, skin irritation, or an unusual reaction occurs, discontinue use immediately and consult your medical professional.  
 Care: Hand wash using mild soap. Rinse thoroughly. Air dry only. Do not tumble dry.

**Air Cell Adjustment:** Air-cells come pre-inflated to correct pressure and do not require adjustment. Change only at high altitudes, or if the normal adjustments are not sufficient to achieve a comfortable fit. To deflate air-cells unfold valves at top of brace. Insert inflation tube 1 1/2" into each air-cell. Air-cell will deflate. To inflate, blow through tube into each air-cell. After adjusting inflation level squeeze valve flat and fold back under the air-cell


**Fitting:** Ankle Stirrup should be worn over an absorbent sock with a shoe. If it does not fit well readjust straps and/or heel pad width

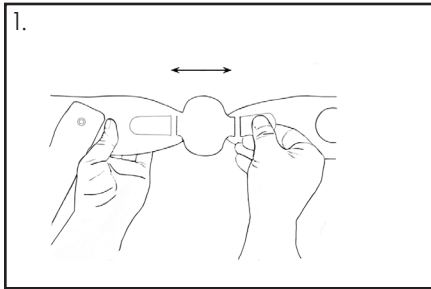


1. Adjust both straps on heel-pad to accommodate heel width.

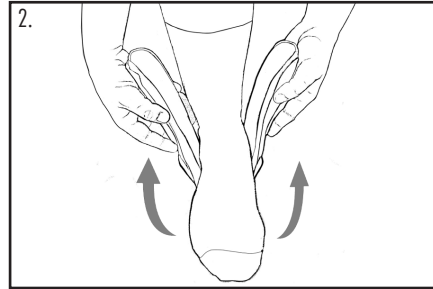


2. Place foot in brace with heel directly on heel-pad and bring the sides up.

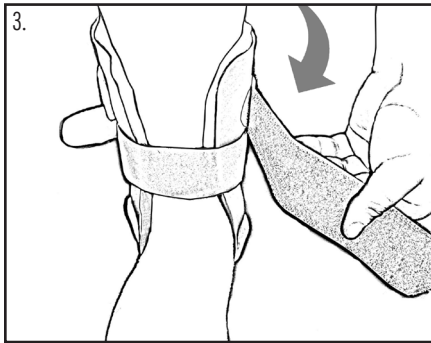

3. Wrap the lower strap first, then wrap the top strap.

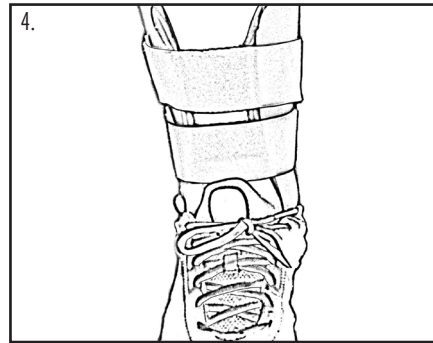

4. Apply shoe over brace.

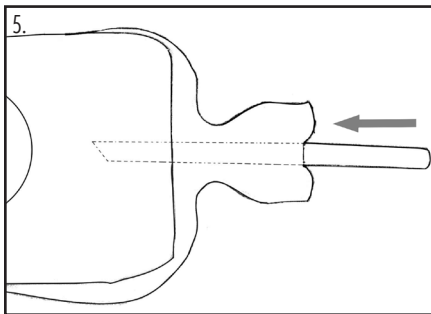

5. Adjust inflation level (if needed)



© 2014 Weber Orthopedics Inc.  
 DBA Hely & Weber all rights reserved  
 1185 E Main St., Santa Paula, California 93060  
 California: 800-221-5465, National: 800-654-3241  
 International: 805-525-4244, fax: 805-933-2348  
 U.S. fax: 800-559-5975, www.hely-weber.com




MDSS GmbH, Schiffgraben 41  
 30175 Hannover, Germany

**Warnings and Instructions: Review carefully, proper application is required**


**⚠ Warning:** This device will not prevent or eliminate risk of injury. Do Not Overtighten. If swelling, pain, skin irritation, or an unusual reaction occurs, discontinue use immediately and consult your medical professional.  
 Care: Hand wash using mild soap. Rinse thoroughly. Air dry only. Do not tumble dry.

**Air Cell Adjustment:** Air-cells come pre-inflated to correct pressure and do not require adjustment. Change only at high altitudes, or if the normal adjustments are not sufficient to achieve a comfortable fit. To deflate air-cells unfold valves at top of brace. Insert inflation tube 1 1/2" into each air-cell. Air-cell will deflate. To inflate, blow through tube into each air-cell. After adjusting inflation level squeeze valve flat and fold back under the air-cell


**Fitting:** Ankle Stirrup should be worn over an absorbent sock with a shoe. If it does not fit well readjust straps and/or heel pad width

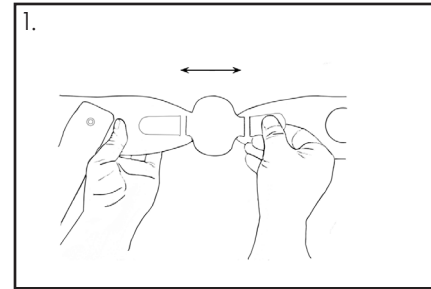


1. Adjust both straps on heel-pad to accommodate heel width.

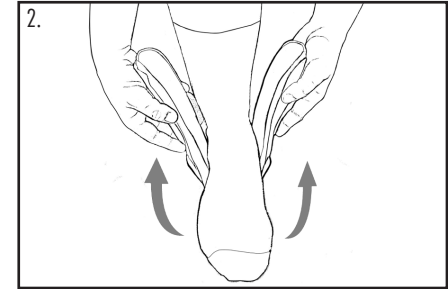


2. Place foot in brace with heel directly on heel-pad and bring the sides up.

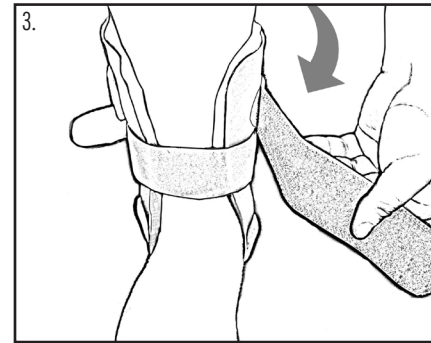

3. Wrap the lower strap first, then wrap the top strap.

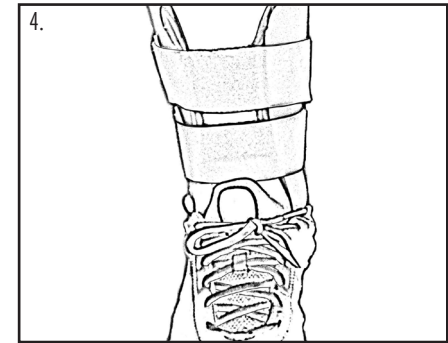

4. Apply shoe over brace.

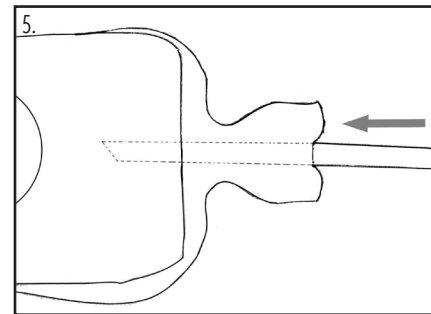

5. Adjust inflation level (if needed)



© 2014 Weber Orthopedics Inc.  
 DBA Hely & Weber all rights reserved  
 1185 E Main St., Santa Paula, California 93060  
 California: 800-221-5465, National: 800-654-3241  
 International: 805-525-4244, fax: 805-933-2348  
 U.S. fax: 800-559-5975, www.hely-weber.com



MDSS GmbH, Schiffgraben 41  
 30175 Hannover, Germany